

Castor

Castor eConsent Site User Guide

Version 2023.3

3
4
6
9
10
12
12
12
16
20
23
26
27
27
29
30
33




1. Register an account

Before you are able to access your eConsent study you will first need to be invited and create an account for eConsent.

Invited users will receive an email with a link.

You've received a new Castor eConsent invitation

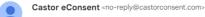

O castor eConsent

You've received a new Castor eConsent invitation

You have been invited to join the study team for 'Castor eConsent' in Castor eConsent. To join the study please click the button below.



By clicking the 'Join the Study' link provided in the email, you will be directed to the registration page:

	S English
	G Light
O Castor eConsent	
Welcome to Castor eConsent	
To get started, please sign up	
First name Last name	
Email	
Create password	
Please use at least:	
× 1 uppercase character	
× 1 numeric character × 1 lowercase character	
× 8 characters	
The second s	
I declare I have read the Terms Of Use	
Sign up	



You will need to accept the Terms of Use and provide the following details when registering the account:

- First name
- Last name
- Email will be automatically pre-filled
- Create password

After filling out the details and accepting the Terms of Use, click on the 'Sign up' button to complete the registration.

Should you already have an account registered in eConsent, you will be directed to the eConsent login page where you can log in with your eConsent credentials.

2. Log In

In order to access your eConsent study after registration has concluded you will need to access the following url's bases on the location of your study, US or Europe:

- US eConsent <u>https://us.castorconsent.com/</u>
- EU eConsent https://eu.castorconsent.com/

Study data is only stored on one of these servers and only one server location can be accessed at a time. If you don't see any studies listed please make sure to check our other server location.

You can toggle between the two study locations by selecting the Europe or USA icons at the top of the page.

: astor eConsent		🙁 🚺 Europe
Search	Q 0	
• Live	The Impact of Yoga and Meditation on Chronic Pain: A Randomize 🛈	Main contact:
• Testing (Castor)	Castor eConsent Study (Testing)	Main contact:
• Testing (Customer)	Castor eConsent Study	
Testing (Customer)	Support Team Study ()	Main contact:



3. Open a study

After logging in, the 'Studies' overview will be displayed with two sections, 'Invitations' and 'Active'.

The 'Invitations' section will only appear if you have been sent an invitation to join a study and have not yet accepted or rejected the invitation.

Before accessing any studies you will first need to accept the invitation by pressing the 'Accept' button.

Rejecting the invitation cannot be undone and you will need to contact the study admin to request a new invite.

S castor eConsent	8 0
	🚺 Europe 📒 USA 🛈
Studies	
Search Q 3	
Invitations	
Testing (Customer) Support Team Study	Accept Reject
Active	
Testing (Castor) Castor eConsent Study	Main contact:
Testing (Customer) Castor eConsent Study	Main contact:
	A A A A A A A A A A A A A A A A A A A
	Help us improve

The 'Studies' overview consists of the following elements:

• Servers: Here the data storage locations are displayed. There are two server locations available: Europe and the United States. Study data is only stored on one of our servers and only one server location can be accessed at a time. If you don't see any studies please make sure to check a different server location.



- Search bar: it is possible to search by organization name, study contact email address or study name
- Invitations: If you are invited to a study, but you have not accepted the invite yet, the study will be listed in the top panel. It is possible to either accept or reject an invitation to a study. Once you have accepted the invitation, you will be able to open a study.
- Active: once you have accepted an invitation to a study, it will appear in the 'Active' list. The list displays all studies you have access to. For each study, the following information is shown:
 - **Status**: Testing (Castor), Testing (Customer), Live. Read more about what each status means in the article: <u>Study status in eConsent</u>
 - Study title: name of the study
 - Main contact: main contact of the study. If the main contact has been invited to the study, but have not accepted the invitation yet, the status will be 'Main contact invited'
 - **Study information**: by clicking on the(i) icon, it is possible to view additional study information

S castor eConsent		8 O
Studies		Europe SUSA ()
Castor eConsent Study	Study Information × Support Team Study	Main contact:
Testing (Customer) Castor eConsent Study Testing (Customer) Support Team Study ①	Study contact details Street: 54 South Crescent City: London Phone number: 079 7693 3810 Email: studyteam@test.com Country: United Kingdom	Main contact:
	Close	Heb us improve



Studies for which an invitation has been accepted will appear in the 'Active' section and it is possible to access the study by clicking on the study row.



4. Adding a participant

To create a new participant in your study, follow the steps below:

- 1. Navigate to the 'Participants' overview
- 2. Press the 'Add participant' button

S castor eConsent			© ©						
Study settings Org & sites Kuser & roles	Participants Search for a participan Search	t © Q						T Filter Add	d participant
	Screening ID $\uparrow\downarrow$	Participant ID 1	Invite status $\uparrow\downarrow$	Eligibility status $\uparrow\downarrow$	Consent status $\uparrow\downarrow$	Invited on $\uparrow\downarrow$	Consented on $\uparrow\downarrow$	ICF version	Site ↑↓
Consent forms	AM50001	•	Active Act	Pending	Awaiting consent	05 Sep 2023 11:24	-	1.0, 01 Aug 2023	Amsterdam Cli
Participants									
③ Audit trail									

- 3. On the 'Add participant' screen, fill out the details:
 - Site: choose a site where a participant should be created.
 - Language: specify the language
 - Add consent form: choose the applicable consent form
 - Send an invitation to participant?: set this option to 'Yes' if you wish to send an email invitation to the participant and to 'No' if no notification emails should be sent.
 - **Participant email address**: if you wish to send an invitation to a participant, make sure to add the participant's email address. If no notification should be sent, please skip this field.

G castor eConsent		The Impact of Yoga and Meditation	on Study State
Study settings	Add participant		
🔠 Org & sites	Screening ID The Screening ID will be generated automatically by eConsent once the participant is created.		
in Har & roles ☆	Site		
Consent forms	Amsterdam Clinic ~		
_	Language		
E Participants	Dutch ~		
3 Audit trail	Participant email address		
	e.g. jess@email.com		
	Send an invitation to participant?		
	Tes No		
	Assign documents		
	These are all the documents available with your selected options. F participant.	Please select the document(s) you want to assig	n to this
	Name	Version	LAR signature
	Consent form - Amsterdam Clinic (Dutch)	1.0, 01 Aug 2023	No
	Add participant Cancel		

If you only have rights for one site or one consent form, the applicable fields will be automatically populated.

- 4. Once you are ready, click on the 'Add participant' button to add the participant. Clicking on 'Cancel' will discard all the information. You will be redirected to the 'participants' overview screen.
- 5. When a participant is created, a screening ID will be generated automatically.

4.1 (Un)Archiving a participant

1. To archive a participant, start by navigating to the 'Participants' listing.

2. Locate the participant you wish to archive, then click on the three-dot context menu in the Options column.

3. From the dropdown menu, select the 'Archive Participant' option. This will securely archive the participant, and please note that after archiving, the participant ID cannot be reused.

castor eConsent				The Impact of	Yoga and Meditation or	n Study Status: Liv	e			®
Study settings	Participants									
Org & sites	Search for a participal	nt ©							∀ Filter	Add participant
h User & roles	Screening ID ↑↓	Participant ID ↑↓	Invite status ↑↓	Eligibility status ↑↓	Consent status ↑↓	Invited on ↑↓	Consented on ↑↓	ICF version	Site ↑↓	Options
	AM50007		Not invited	Pending	-	-		1.0, 01 Aug 2023	Amsterdam Clinic	
Consent forms	AMS0006		🖸 Invited	Pass	Awaiting consent	17 Sep 2023 11:44		1.0, 01 Aug 2023	Amsterdam Clinic	:
Participants	AM50005		🖸 Invited	Pending	Awaiting consent	06 Sep 2023 18:04		1.0, 01 Aug 2023	Amsterdam Clinic	:
p Participants	BER0001		Not invited	Pass	Consented		01 Sep 2023	1.0, 09 Aug 2023	Berlin Clinic	:
) Audit trail	AM50004	-	Not invited	Pending	Awaiting consent			1.0, 01 Aug 2023	Amsterdam Clinic	1
Addit trail	AM50003		Not invited	Pass	Consented		05 Sep 2023	1.0, 01 Aug 2023	Amsterdam Clinic	:
	AMS0002	-	Not invited	Pass	-			1.0, 01 Aug 2023	Amsterdam Clinic	:
	AMS0001	-	Ø Active	Pass	Awaiting consent	05 Sep 2023 11:24	-	1.0, 01 Aug 2023	Amsterdam Clinic	:
									A	rchive

4. Once the participant ID is archived, it will no longer be displayed in the 'Participants' page. To locate the archived participant ID, click on the 'Filter' button and choose 'Include archived' option, afterwards click 'Apply filters':

astor eConsent			The Impact of	Yoga and Meditat	ion on Study Status: L	ive			8
	AMS0003	-	Not invited	Pass	Consented		05 Se	End	Ħ
Study settings	AM50002	-	Not invited	Pass	-	-		Eligibility status	
Org & sites								All × ×	•
User & roles								Consent status	
								All × ×	•
Consent forms								Consented betwee	een
								Start	Ħ
Participants								End	Ħ
A								Countries	
Audit trail									•
								Site	
								All × ×	•
								Include arc	nived
								Apply filter	;
								Clear all	

5. Locate the archive participant, click on the context menu (three dots) and choose 'Unarchive' option:

Castor eConsent			The Ir	npact of Yoga and Med	itation on Study S	itatus: Live			8
Study settings	Participan	ts							
Org & sites	Search for a particular Search		Q			▼ Filter A	Add participant		
¦ி User & roles	ility status ↑↓	Consent status ↑↓	Invited on ↑↓	Consented on ↑↓	ICF version	Site ↑↓	Options		
Consent forms	nding	Consented	17 Sep 2023 14:12	04 Oct 2023 10:39	1.0, 01 Aug 2023	Amsterdam Clinic	:	∀ Filter >	<
Consent forms	nding	Consented	-	10 Oct 2023 17:19	1.0, 01 Aug 2023	Amsterdam Clinic	:	Invite status	
Participants	s	Awaiting consent	17 Sep 2023 11:44		1.0, 01 Aug 2023	Amsterdam Clinic	:	All × × v	
	nding	Awaiting consent	06 Sep 2023 18:04	-	1.0, 01 Aug 2023	Amsterdam Clinic Berlin Clinic	:	Invited between	
S Audit trail	s	Consented Awaiting consent	-	01 Sep 2023	1.0, 09 Aug 2023	Amsterdam Clinic	:	Start	3
	s	Consented		05 Sep 2023	1.0, 01 Aug 2023	Amsterdam Clinic	:	End	E
	s				1.0, 01 Aug 2023	Amsterdam Clinic	:	Eligibility status	
	s	Awaiting consent	05 Sep 2023 11:24	-	1.0, 01 Aug 2023	Amsterdam Clinic	:	All × × ×	
						Ur	narchive	Consent status	
								All × × ×	
								Consented between	
								Start É	E E



5. Signing a eConsent form

Castor eConsent platform supports two consent methods: Electronic (created with the eConsent ICF Builder) or Paper (used to upload a form signed on paper). Consenting can be done electronically with a participant account or without a participant account (in-person signing). In cases where a paper consent method is used, signed consent forms can be uploaded into the eConsent platform.

If your desired consent form is not visible or your desired consent method (paper or electronic) is not present, please reach out to the study admin and request the study admin to add the consent method or consent form to your site.

5.1 Electronic signing

When choosing an electronic consent method, participants can consent online with an account or in-person without creating an account.

5.1.1 Electronic signing with an account (remotely)

To initiate the electronic signing, follow these steps:

- 1. Add a participant from the 'Participants' tab
- 2. When adding a new participant, choose 'Electronic' consent method
- 3. In the 'Participant email address' add the participant's email
- 4. If the signing procedure happens remotely, select 'Yes' in the 'Send an invitation to participant?' field
- 5. Choose the relevant documents in the 'Assign documents' section
- 6. Click on the 'Add participant' button create a participant and send an invite



S castor eConsent		The Impact of Yoga a	nd Meditation on	Study Status: Live
Study settings	Add participant			
🔠 Org & sites	Screening ID			
ന്ന് User & roles	The Screening ID will be generated automatically by eConsent once the participant is created.			
	Site			
Consent forms	Amsterdam Clinic ~			
	Language			
Participants	Dutch ~			
🕚 Audit trail	Consent method Electronic Paper (used to upload a form signed on paper)			
	Participant email address			
	e.g. jess@email.com			
	Send an invitation to participant?			
	Ves No			
	Assign documents			
	These are all the documents available with your selected options. Ple participant.	ease select the document(s) you want to assig	gn to this	
	Name	Version	LAR signature	
	Consent form - Amsterdam Clinic (Dutch)	1.0, 01 Aug 2023	No	
	Add participant Cancel			


- 7. You will then be redirected to the 'Participants' listing.
- 8. From the 'Participants' listing, double click on the participant row
- 9. You will be presented with the participant details
- 10. If eligibility is checked in the study, the 'Eligibility status' will be automatically set to 'Pending'. Click on the 'Update eligibility status' button to change the eligibility status. Choose from the following options: Pass, Pending, Fail. If an integration is set up with a Castor CDMS study, the eligibility field might be updated automatically, check with the study admin if this is applicable for your study.
- 11. To allow a participant to sign the form, in the 'Consent information (most recent)' section, click on the 'Unlock' button. Please confirm you would like to update the form status to unlocked in the pop-up window. The form will now become available for the participant for signature.
- 12. In order to prompt your invitee once more, simply click on the 'Send Reminder' button to resend the invitation.



A participant can electronically sign a consent form by creating an account from the invitation email. When a participant clicks on the 'Join the study' button, they need to <u>create an account</u> if they don't have one yet. Participants with an account, can login using their credentials.

٢	You've received a new Castor eConsent invitation								
	S castor eConsent								
	Netherlands Institute has invited you to their study Netherlands institute has invited you to participate in a research study named "castor excesser. Before joining his research study you are asked to sign a consent form. To review the study and consent form, please click the button below. Join the study								
	us.castorconsent.com Faster, smarter medical research								

After setting up their account and logging in, a dashboard will be presented.

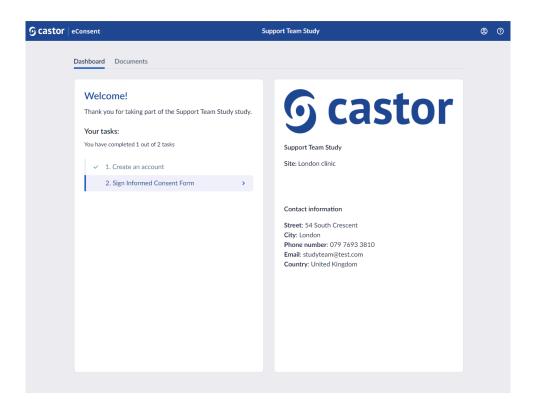



The participant will need to Press the 'Sign Informed Consent Form' where they can download the consent form or view and sign.

ତ castor	eConsent	Support Team Study						
43	Dashboard	Documents						
	Informed Consent Forms							
		Consent Form London clinic ICF 1.0, 01 Aug 2022 Awaiting consent View & Sign						

Upon selecting 'View & Sign', the participant will then select the required fields, type their full name, add their electronic signature, then press the 'Submit' button to submit their signature.

€ 0	onsent form: London clinic ICF 1.0, 01 Aug 2022	C Download
	All information taken from the study will be coded to protect each subject's name. No names or other identifying information will be used when discussing or reporting data. The investigator(s) will safely keep all files and data collected in a secured locked cabinet in the principal investigators office. Once the data has been fully analyzed it will be destroyed.	
	"Required indicates required checkbox I authorize the use of my records, any observations, and findings found during the course of this study for education, publication and/or presentation "Required I voluntarily agree to participate in this research program "Required Subject signature	
	Date 06 October 2022 Full name e.g. Jess Smith Signature	
	Submit Clear	

5.1.2 Electronic signing via a Video-call (remotely)

When the option 'Video calling for signing' is enabled in the study settings tab, it is possible to obtain informed consent remotely via a video call.


1. Once the configuration is enabled, when opening an individual participant, a 'Consent call' card will appear:

G castor eConsent			The Impact of Yoga and Me	editation on Study Status: Live
Study settings	Profile Consent Forms	Audit trail		
🔠 Org & sites	Participant profile			Personal details
🏦 User & roles	Screening ID	Participant ID	Created by	<u> </u>
Consent forms	AMS0001 Site	- Language	Tonya Support Invite status	Date of birth Gender
Participants	Amsterdam Clinic	Dutch	② Active	Contact information
③ Audit trail	Consent information (most	recent)	Send reminder	Consent call
	Consent status Awaiting consent	Consent form sent on 05 Sep 2023 11:24	Consent form signed on	Start call
	Consent form	Version	Consent method	
	Consent form - Amsterda m Clinic (Dutch)	1.0, 01 Aug 2023	Electronic	
	Complete information			
	Screening information		Update eligibility status	
	Eligibility status Pass	Changed on 05 Sep 2023 14:52		

 Click on the 'Start call' button to initiate the video call. A pop-up window will appear which allows to specify whether an ICF form should be included in the call. Specify which ICF should be included in the 'Select ICF' dropdown field and click 'Start call' button to proceed to the call:

G castor eConsent	nt The Impact of Yoga and Meditation on Study Status: Live						
Study settings	← Profile Consent Forms	Audit trail					
🔠 Org & sites	Participant profile				Personal details		
🏦 User & roles	Screening ID	Participant ID	Created by		© -		
Consent forms	AMS0001 Site	- Language	Tonya Support Invite status		Date of birth Gender		
Participants	Amsterdam Clinic	Dutch	③ Active		Contact information		
③ Audit trail	Consent information (most	recentl	video call		*		
	Consent status	Consent fc	Name	Versie	on Consent status		
	Awaiting consent	05 Sep 20	Consent form - Amsterdam Clinic (Dutch) 1.0, (01 Aug 2023 Awaiting consent		
	Consent form Consent form - Amsterda m Clinic (Dutch)	Version 1.0, 01 Au			Cancel Start call		
	Screening information		Update eligibility sta	tus			
	Eligibility status Pass	Changed on 05 Sep 2023 14:52					

- 3. You will be redirected to the call screen. During the call, the following options are presented:
 - Unlock to enable signing/Lock to disable signing which allows to enable the ICF for signing and lock it after the form has been signed
 - Start video/Stop video which allows to turn on/off the investigator's camera
 - Mute/Unmute to turn on/off the investigator's microphone
 - End call button to finish the call
 - Arrow to return back to the participant overview

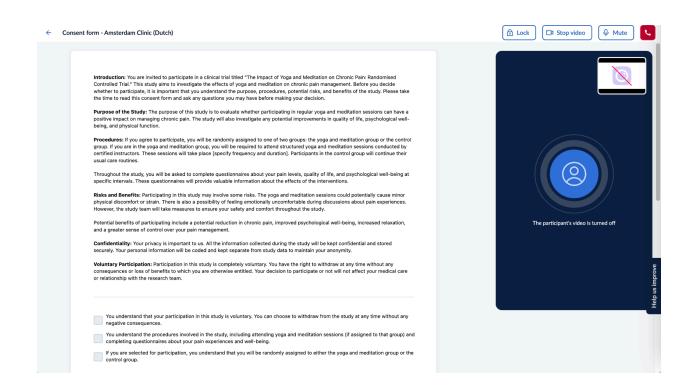

4. The participant will need to login to the participant portal and choose the relevant study in the 'Studies' overview.

5. Once the investigator starts the call, the participant will see an invitation to 'Join call' within the 'Dashboard' tab.

Scastor eConsent	Тһ	e Impact of Yoga and Meditation on	0	?
	Dashboard Documents Study Profile			
	Your call has started, please join the call.	Cit Join call		
	Welcome! Thank you for taking an interest in this study. Your tasks: You have completed 1 out of 2 tasks ✓ 1. Create an account 2. Sign Informed Consent Form Scheduled video call You can join the call The study team member started the call, please use the button below to join the call. Image: Join call	The Impact of Yoga and Meditation on Chronic Pain: A Randomized Controlled Trial Site: Amsterdam Clinic Chronic pain is a complex and pervasive health issue affecting millions of individuals worldwide. Traditional medical interventions often have limitations in providing heisike relief, which has let on a increased interest in complementary and alternative therapies. This randomized controlled Read more Contact information Medfech Innovate 456 Health Drive, Suite 789 Medfech Texability Bit 23 CD United Kingdom		Help us improve

5. To finish a call, click on the End call button



Please refer to the article <u>Troubleshooting Video/Audio</u> on how to resolve potential video/audio issues.



5.1.3 Electronic signing in-person

If a participant is invited to sign the consent form online but instead wishes to sign in-person, if the study allows, follow the steps below.

In-person signing, with or without an account

If the study allows in-person signing, the study team members with the 'Site Admin', 'Study Investigator', 'Site Investigator' roles can start the in-person signing flow:

- 1. Navigate to the 'Participants' overview and click on the 'Add participant' button or choose an existing participant.
- 2. When creating a new participant, you may want to set the field 'Send an invitation to participant' to 'No'. In this case, an email invite will not be sent and the invite status in the Participants listing will be set to 'Not invited'.

S castor eConsent		The	Impact of Yoga a	nd Meditation on
	Add participant			
Study settings	Screening ID			
🔠 Org & sites	The Screening ID will be generated automatically by eConsent once the participant is created.			
ີ່ User & roles	Site			
Consent forms	Amsterdam Clinic ~			
	Language			
Participants	Dutch ~			
	Consent method			
S Audit trail	O Electronic			
	Paper (used to upload a form signed on paper)			
	Participant email address			
	e.g. jess@email.com			
	Send an invitation to participant?			
	Yes 🔾 No			
	Assign documents			
	These are all the documents available with your selected options participant.	s. Please select the document(s) you want to assig	gn to this
	Name	Ver	rsion	LAR signature
	Consent form - Amsterdam Clinic (Dutch)	1.0), 01 Aug 2023	No
	Add participant Cancel			




3. Participants who received an invite per email can also use the in-person signing flow. In this case, the option to sign in-person will appear in the 'Consent information (most recent)' section in the context menu (three dots menu):

G castor eConsent			The Impact of Yoga and Med	litation on Study Status: Live	© ©
Study settings	Profile Consent Forms	Audit trail			
🖾 Org & sites	Participant profile			Personal details	
n User & roles	Screening ID AM50006	Participant ID	Created by	® -	
Consent forms	Site	- Language	Tonya Support Invite status	Date of birth Gender 	
Participants	Amsterdam Clinic	Dutch	☑ Invited	Contact information	
🕚 Audit trail	Consent information (most	recent)	Send reminder :	Consent call	
	Consent status Awaiting consent	Consent form sent on 17 Sep 2023 11:44	View Conser - Download	Start call	
	Consent form Consent form - Amsterda m Clinic (Dutch)	Version 1.0, 01 Aug 2023	Conser Lock Electro Sign in-person Add new		
	Screening information		Update eligibility status		Help us improve
	Eligibility status Pass	Changed on 17 Sep 2023 11:47			Hel

4. In case where no invite to a participant is sent and when there is already an ICF added which hasn't been signed by the participant yet (status of the ICF is "Awaiting consent" if the invite is already sent, or "undefined" if there is no invite sent), in the 'Consent information (most recent)' section, click on the "Sign in-person" option:

GCASTOR eConsent			The Impact of Yoga and Meditation on	- Study Status: Live	
Study settings	← Profile Consent Forms	Audit trail			
Org & sites	Participant profile			Personal details	
n User & roles	Screening ID	Participant ID	Created by	Ø -	
Consent forms	AM50002	-	Tonya Support	Date of birth	Gender
	Site	Language	Invite status		-
Participants	Amsterdam Clinic	Dutch	Not invited	Contact information	
Audit trail	Consent information (most	recent)	Sign in-person		
				Consent call	
	Consent status	Consent form sent on	Consent form signed on		
		-		Start call	
	Consent form	Version	Consent method		
	Consent form - Amsterda m Clinic (Dutch)	1.0, 01 Aug 2023	Electronic		

5. After clicking on the "Sign in-person" options, a confirmation dialog will be shown with the following message: "The form will be unlocked for signing, the study team member will be required to monitor the consent form so that the signer does not access other parts of the application." Select the checkmark for the confirmation statement "I have verified the signer's identity as part of the in-person signing procedure." and click on the 'Confirm identification' button to proceed further:

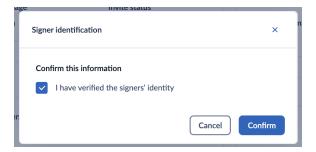



5.2 Consenting on paper

If the ICF template uses a paper consent method and this method is enabled in the study Settings tab, a scan/copy of the paper signed ICF can be uploaded against the template. Make sure that a paper consent form is created in the 'Consent forms' tab first, otherwise it will not be possible to upload a paper ICF.

1. When adding a participant, choose the 'Paper (used to upload a form signed on paper)' and select the applicable paper consent form.

ら castor eConsent	The Impact of Yoga and Meditation on Study Stat
~	
Study settings	
Add participant	
Muser & roles Screening ID The Screening ID will be generated automatically by eConsent once the participant is created.	
Consent forms Site	
Amsterdam Clinic ~	
Participants Language	
Dutch	
Audit trail Consent method	
Electronic	
Paper (used to upload a form signed on paper)	
Assign documents	
These are all the documents available with your selected optio participant.	ns. Please select the document(s) you want to assign to this
Name	Version
Consent form - Amsterdam clinic (Paper - Dutch)	1.0, 01 Aug 2023
Add participant Cancel	

- 2. After the participant is added, in the participants overview, click on the participant row.
- In the 'Consent information (most recent)' section, click on the 'Upload signed form' button.

castor eConsent			The Impact of Yoga and Meditation on	- Study Status: Live	
Study settings	Profile Consent Forms	Audit trail			
Org & sites	Participant profile			Personal details	
User & roles	Screening ID	Participant ID	Created by	® -	
Consent forms	AM50003 Site	- Language	Tonya Support Invite status	Date of birth -	Gender -
Participants	Amsterdam Clinic	Dutch	😪 Not invited	Contact information	
Audit trail	Consent information (most	t recent)	Upload signed form		
	Consent status	Consent form sent on	Consent form signed on		
	Awaiting consent	-	-		
	Consent form	Version	Consent method		
	Consent form - Amsterdam clinic (Paper - Dutch)	1.0, 01 Aug 2023	Paper		

4. In the 'Upload signed form' pop-up window, confirm the requested details and upload the signed form. To save the form, click on the 'Confirm' button.

G castor eConsent			The Impact of Yoga and Meditation on Study Status: Live			© ©
Study settings	← Profile Consent Forms	Audit trai	I			
🔠 Org & sites	Participant profile	16	Upload signed form	ails		
in the second se	Screening ID	Partici		-		
Consent forms	AMS0003 Site	- Langua	Confirm this information I have verified the signers' identity		Gender -	
Participants	Amsterdam Clinic	Dutch	This document has been signed by all parties	nation		
① Audit trail	Consent information (most	t recent)	Signed on 05.09.2023			
	Consent status Awaiting consent	Consen	Upload your signed form below ⊘ signed_icf.pdf □			
	Consent form Consent form - Amsterdam clinic (Paper - Dutch)	Version 1.0, 01	 I confirm this is Consent form - Amsterdam clinic (Paper - Dutch) 1.0, 01 Aug 2023 			Help us improve
			Cancel			Help us
	Screening information		Update eligibility status			
	Eligibility status Pass	Changed 05 Sep 20				

Accepted files: .pdf, .jpeg, .jpg, .png and the maximum file size is 15 Mb.



5. After the file is uploaded, the 'Consent status' will be automatically changed to 'Consented' in the 'Consent information (most recent)' section in case there are no other pending documents for signature.

G castor eConsent		т	he Impact of Yoga and Meditation on	- Study Status: Live	© (2)
Study settings	← Profile Consent Forms	Audit trail			
🔠 Org & sites	Participant profile			Personal details	
📸 User & roles	Screening ID	Participant ID	Created by	® -	
Consent forms	AMS0003 Site	- Language	Tonya Support Invite status	Date of birth Gender 	
E Participants	Amsterdam Clinic	Dutch	😪 Not invited	Contact information ☑ -	
🕚 Audit trail	Consent information (most	recent)	Add new :		
	Consent status	Consent form sent on	Consent form signed on		
	Consented	-	05 Sep 2023		
	Consent form Consent form - Amsterda m clinic (Paper - Dutch)	Version 1.0, 01 Aug 2023	Consent method Paper		Help us improve
	Screening information		Update eligibility status	-	He

6. If you would like to view, download or remove the uploaded file, click on the options menu (three dots) in the 'Consent information (most recent)' section:

Scastor eConsent			The Impact of Yoga and Meditation on	- Study Status: Live	
Study settings	← Profile Consent Forms	Audit trail			
🔠 Org & sites	Participant profile			Personal details	
💏 User & roles	Screening ID	Participant ID	Created by	© -	
Consent forms	AMS0003 Site	- Language	Tonya Support Invite status	Date of birth	Gender -
E Participants	Amsterdam Clinic	Dutch	🔓 Not invited	Contact information M -	
🕄 Audit trail	Consent information (most	recent)	Add new 🚺		
	Consent status Consented	Consent form sent on	View Consen 05 Sep Download		
	Consent form	Version	Consen Remove upload		
	Consent form - Amsterda m clinic (Paper - Dutch)	1.0, 01 Aug 2023	Paper		



If a participant/LAR has been invited to the study, they will also be able to view the uploaded paper form.

5.3 Print-to-sign

The ICF builder for ICF templates that use a paper consent method allows pdf file uploads. This means that the finalized "paper" ICF file can be uploaded against the eConsent paper ICF template. By doing so, the actual content of the ICF becomes available for site members. The site member can download and print the ICF template directly from the eConsent participant record, in order to make sure the right ICF version is getting signed.

1. To download, open the participant record and in the "Consent Information (most recent) dialogue box select the options menu (three dots) to the right and select download.

castor eConsent		т	he Impact of Yoga an	d Meditation on •	- Study Status: Live	
 Study settings Org & sites 	Profile Consent Forms Participant profile	Audit trail			Personal details	
User & roles	Screening ID BER0001	Participant ID -	Created by Tonya Support Invite status		③ - Date of birth	Gender
Participants	Site Berlin Clinic	Language English	Invite status	1	Contact information	
udit trail	Consent information (most	recent)	Upload signed			
	Consent status Awaiting consent	Consent form sent on	Consent forr -	View		
	Consent form Consent Form - Berlin Cli nic (English - Paper)	Version 1.0, 09 Aug 2023	Consent met Paper	Add new		

2. Once the printed version is signed by all parties, a scan can be uploaded into the system by selecting "Upload signed form".



3. In the dialog box confirm the required information and upload the file by selecting "choose file" or drag and drop the file into the box and confirm.

Accepted file formats: .pdf, .jpeg, .jpg, .png and the maximum file size is 15 Mb.

5.4 Reconsenting

If a participant needs to be reconsented during the study, follow the steps below:

- 1. Select the participant row to access the participant's profile.
- 2. Press the 'Add New' button to assign the new consent form.

Castor eConsent		т	he Impact of Yoga and Meditation on	- Study Status: Live	
Study settings	← Profile Consent Forms	Audit trail			
Org & sites	Participant profile			Personal details	
h User & roles	Screening ID	Participant ID	Created by	@ -	
	BER0001	-	Tonya Support	Date of birth	Gender
Consent forms	Site	Language	Invite status	-	-
	Berlin Clinic	English	😪 Not invited	Contact information	
Participants				⊠ -	
Audit trail	Consent information (most	recent)	Add new :		
	Consent status	Consent form sent on	Consent form signed on		
	Consented	-	01 Sep 2023		
	Consent form	Version	Consent method		
	Consent Form - Berlin Cli nic (English - Paper)	1.0, 09 Aug 2023	Paper		

- 3. In a new pop-up screen, you can repeat the process for on-site, remote or paper consent by selecting the relevant consent method and consent form.
- 4. Once required settings are set and the consent form is selected, press 'Add' to continue.
- 5. The consent status will be updated to 'Awaiting consent'.

6. Exporting consent form and participant information

Users with the 'Study Investigator' or 'Study Monitor' role are able to generate the export that contains consent form signature or participant personal details information.



To export the consent form or personal details information, follow the steps below:

G castor eConsent			Castor eConsent Study	- Study Status: Testing	(Customer)		© (2
Study settings	Subjects						_
🔠 Org & sites	Search for a subject ③					Exp	ort 7 Filter
	Search	Q					
% User & roles							
	Screening ID ↑↓	Subject ID ↑↓	Invite status $\uparrow\downarrow$	Eligibility status $~\uparrow\downarrow~$	Consent status $\uparrow\downarrow$	Invited on $\uparrow\downarrow$	Consented on $\uparrow\downarrow$
Consent forms	TEST-LCL0006	-	Active Act	Pass	Consented	09 Nov 2022 12:00	14 Nov 2022 13:27
	TEST-LCL0005	-	Active Act	Pass	Awaiting consent	09 Nov 2022 11:26	
E: Subjects	TEST-LCL0004	-	🖾 Invited	Pass	Awaiting consent	09 Nov 2022 10:52	
	TEST-LCL0003	-	Not invited	Pending	Awaiting consent	-	
S Audit trail	TEST-LCL0002	-	😪 Not invited	Pending	-	-	· ·
	TEST-LCL0001	-	😪 Not invited	Pending	Awaiting consent	-	


1. In the 'participants Overview', click on the 'Export' button:

2. In the 'Generate export file' pop-up window, in the 'Export type', select 'Consent form information' or 'Participant personal details' option. Click on the 'Generate' button to proceed with the export. Only signature information for which you have view rights will be included in the export file.

Study settings	Participants					
Org & sites	Search for a participa	nt 0			Exp	ort 7 Filter
user & roles	Screening ID ↑↓	Participant ID	Generate export file ×	Invited on ↑↓	Consented on ↑↓	ICF version
Consent forms	AM50008	-	Export type	17 Sep 2023 14:12	04 Oct 2023 10:39	1.0, 01 Aug 2023
	AM50007	•	Consent form information ^	- 17 Sep 2023 11:44	10 Oct 2023 17:19	1.0, 01 Aug 2023
Participants	AM50005		Consent form information Participant personal details	06 Sep 2023 18:04	•	1.0, 01 Aug 2023
Audit trail	BER0001				01 Sep 2023	1.0, 09 Aug 2023
	AMS0004 AMS0003		Only information for which you have view rights will be included in the export file.		- 05 Sep 2023	1.0, 01 Aug 2023 1.0, 01 Aug 2023
	AMS0002			•	•	1.0, 01 Aug 2023
			Cancel Generate			



3. Once the export is generated, the requester is notified by an email that includes a personal download link. The generated export can then be downloaded from within eConsent by clicking on the 'Download' button.

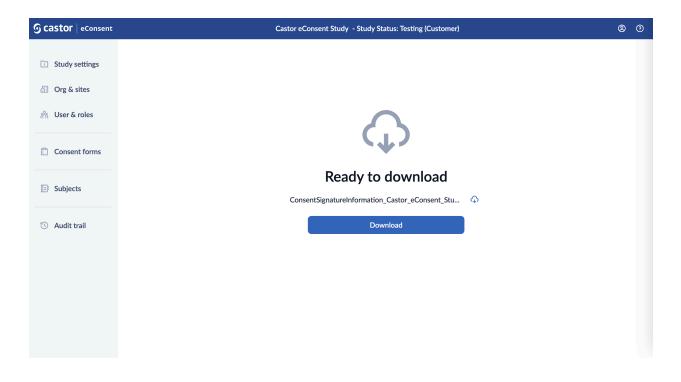


7. Audit trail

The audit trail is a log of study related events:

Castor eConsent			Support Team Study - S	Study Status: Testing (Customer)		8
Study settings	Audit trail						
Org & sites	Date and Time	Actor	Type of user	Screening Id	Action	Additional Info	
	04-Oct-2022 09:47:38 UTC		Study team member	TEST-LND0001	Subject viewed	['siteName':'London clinic',"abbreviation':'LND',"countryCode':'	"GBR"
☆ User & roles	04-Oct-2022 09:47:38 UTC		Study team member	TEST-LND0001	Consent form/information view	w ['siteName':'London clinic',"abbreviation':''LND',"countryCode':'	"GBR"
	04-Oct-2022 09:03:10 UTC		Study team member	TEST-LND0001	Paper signed ICF uploaded	['siteName':'London clinic',"abbreviation':''LND',"countryCode':'	"GBR"
Consent forms	04-Oct-2022 09:03:10 UTC		Study team member	TEST-LND0001	Consent form/information vie	w ['siteName':'London clinic',"abbreviation':''LND'','countryCode':'	"GBR"
	04-Oct-2022 09:01:25 UTC	and the second se	Study team member	TEST-LND0001	Consent form/information view	w ['siteName':'London clinic',"abbreviation':''LND',"countryCode':	"GBR
Subjects	04-Oct-2022 09:01:24 UTC		Study team member	TEST-LND0001	Subject edited	["siteName":"London clinic","abbreviation":"LND","countryCode":"	"GBR
	04-Oct-2022 09:01:24 UTC		Study team member	TEST-LND0001	Subject invited to sign ICF	['siteName':'London clinic',"abbreviation':''LND',"countryCode':'	"GBR
3 Audit trail	04-Oct-2022 09:01:24 UTC		Study team member	TEST-LND0001	Subject viewed	["siteName":"London clinic";"abbreviation":"LND";"countryCode":"	"GBR
	04-Oct-2022 08:52:18 UTC	and the second se	Study team member	TEST-LND0001	Consent form/information vie	w ["siteName":"London clinic";"abbreviation":"LND";"countryCode":"	"GBR
	04-Oct-2022 08:52:17 UTC		Study team member	TEST-LND0001	Subject viewed	['siteName":"London clinic";"abbreviation":"LND";"countryCode":"	"GBR
	04-Oct-2022 08:52:05 UTC		Study team member	TEST-LND0002	Consent form/information view	w ("siteName":"London clinic","abbreviation":"LND","countryCode":"	:"GBR
	04-Oct-2022 08:52:04 UTC		Study team member	TEST-LND0002	Subject viewed	["siteName":"London clinic";"abbreviation":"LND";"countryCode":"	"GBR
	28-Sep-2022 11:58:09 UTC	and the second se	Subject	TEST-LND0007	Consent form/information vie	w ["siteName":"London clinic";"abbreviation":"LND";"countryCode":"	"GBR
	28-Sep-2022 11:57:15 UTC		Subject	TEST-LND0007	Subject accepted invitation	['siteName':*London clinic'**abbreviation':*LND'**countryCode':*	"GBR
	28-Sep-2022 11:54:15 UTC		Study team member	TEST-LND0007	Subject record created	["siteName":"London clinic";"abbreviation":"LND";"countryCode":"	"GBR
	28-Sep-2022 11:54:15 UTC		Study team member	TEST-LND0007	Subject invited to sign ICF	["siteName":"London clinic","abbreviation":"LND","countryCode":"	"GBR
	28-Sep-2022 10:49:23 UTC		Study team member		Study updated	0	
	28-Sep-2022 10:44:13 UTC		Study team member	TEST-LND0001	Subject edited	["siteName":"London clinic";"abbreviation":"LND";"countryCode":"	"GBR
	28-Sep-2022 10:44:13 UTC		Study team member	TEST-LND0001	Subject viewed	["siteName":"London clinic","abbreviation":"LND","countryCode":"	"GBR
	28-Sep-2022 10:44:13 UTC		Study team member	TEST-LND0001	Consent form/information view	w ["siteName":"London clinic","abbreviation":"LND","countryCode":	"GBR
	28-Sep-2022 10:44:03 UTC		Study team member	TEST-LND0001	Consent form/information vie	w ["siteName":"London clinic","abbreviation":"LND","countryCode":"	"GBR
	28-Sep-2022 10:44:02 UTC		Study team member	TEST-LND0001	Subject viewed	["siteName":"London clinic";"abbreviation":"LND";"countryCode":"	"GBR

To view additional details for each event, click on the row:

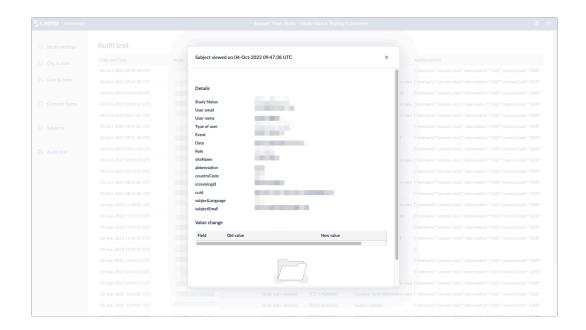

7.1 Audit trail export

Study roles with PII access (i.e. study investigators/monitors) will see an 'Export' button on the global audit trail page they can use to run the export.

Castor eConsent			The Impact of Yoga	and Meditation on S	tudy Status: Live		8
Study settings	Audit trail						Export
Org & sites	Date and Time	Actor		Type of user	Screening Id	Action	Additional Info
User & roles	04-Jan-2024 07:15:25 UTC			Study team member		Study team member edited	{"siteName":"Amsterdam Clinic
User & roles	04-Jan-2024 07:05:12 UTC			Study team member	AMS0006	Consent form/information view	{"siteName":"Amsterdam Clinic
Consent forms	04-Jan-2024 07:05:11 UTC			Study team member	AMS0006	Consent form/information view	{"siteName":"Amsterdam Clinic
	04-Jan-2024 07:05:11 UTC			Study team member	AMS0006	Participant viewed	{"siteName":"Amsterdam Clinic
	04-Jan-2024 07:05:01 UTC			Study team member	AMS0006	Participant viewed	{"siteName":"Amsterdam Clinic
E Participants	04-Jan-2024 07:05:00 UTC			Study team member	AMS0006	Consent form/information view	{"siteName":"Amsterdam Clinic
	18-Dec-2023 12:55:05 UTC			Study team member	AMS0006	In-person consent procedure st	{"siteName":"Amsterdam Clinic
Audit trail	18-Dec-2023 12:55:05 UTC			Study team member	AMS0006	Consent form/information view	{"siteName":"Amsterdam Clinic
	18-Dec-2023 12:55:00 UTC			Study team member	AMS0006	Consent form/information view	{"siteName":"Amsterdam Clinic
	18-Dec-2023 12:54:59 UTC			Study team member	AMS0006	Participant viewed	{"siteName":"Amsterdam Clinic
	13-Dec-2023 08:08:24 UTC			Study team member	AMS0006	Participant viewed	{"siteName":"Amsterdam Clinic
	13-Dec-2023 08:08:14 UTC			Study team member		Study team member accepted	{"email":
	20-Nov-2023 09:40:41 UTC			Study team member		Study team member invited	{"email":
	20-Nov-2023 09:40:26 UTC			Study team member		Study team member removed	{"email":
	20-Nov-2023 09:09:53 UTC			Study team member		Study team member invited	{"email":
	20-Nov-2023 09:09:53 UTC			Study team member		Study team member invited	{"email":
	20-Nov-2023 09:09:53 UTC			Study team member		Study team member invited	{"email":

A pop-up window will appear, click 'Generate' to run the export:

	an at a word.		Church a barran mananah an	10000144
וח	Genera	te export file		×
ny				
ny ny	0	Only audit trail events fo will be included in the ex		rights
Ŋ	Export	format		
ny	-	lionnat		
	Excel			
ny				
ny			Cancel	enerate
O				

Only audit trail events for which you have view rights will be included in the export file.

You will receive an email notification when the file is ready to download. Click the link in



the email or 'Download' button to download the audit trail export. The Export is provided in the Excel format.

G castor eConsent	
Your export is ready for download	
Hello	
Your requested export file of AuditTrailInformation_The_ Impact_of_Yoga_and_Meditation_on_Chronic_Pain:_A_Randomized_ Controlled_Trial_04-01-24_07-19-34.xlsx is now available.	
Click the following link https://eu.castorconsent.com/studies	or use the
button below to download.	-
Download	

You will then be redirected to the download page, click the 'Download' button

Scastor eConsent	The Impact of Yoga and Meditation on Study Status: Live	8	0
Study settings			
Org & sites			
🟦 User & roles	$\mathbf{\circ}$		
Consent forms	(†)		
E Participants	Ready to download		
🕚 Audit trail	AuditTrailInformation_The_Impact_of_Yoga_and_Me		
	Download		
			Help us improve

8. Further Information

If you would like to view our Castor eConsent video tutorials, you can do so here.

For more information regarding creating and managing study as a study admin, please check Castor EDC's knowledge base: <u>https://helpdesk.castoredc.com/en_US/castor-econsent-manual</u>. If you have any questions or concerns, please contact us at <u>support@castoredc.com</u>